

Anti-pNFkBp65 [S529] (K10x)-166Er

Pathologist-Verified Clone for Imaging Mass Cytometry™

Catalog number: 3166026D

Package size and concentration: 25 µg, 0.5 mg/mL

Clone: K10-895.12.50

Isotype: Mouse IgG2b

Pathologist-verified on: Human FFPE

Fluidigm tested on: Human FFPE, Mouse FFPE

Reported reactivity: Human, Mouse

Formulation: Antibody stabilizer with 0.05% sodium azide

Storage: Store at 4 °C. Do not freeze.

Application: IMC paraffin

Technical Information

Description: NFκB-p65, also known as RelA, is a subunit of the NFκB transcription factor complex, an important mediator of inflammatory and immune responses. NFκB-p65 is normally sequestered in the cytoplasm through an interaction with IκB. NFκB translocates to the nucleus when IκB is degraded in response to stimuli such as TNFα. In mammals, 5 members of the NFκB family can form heterodimers in the nucleus to activate different sets of genes. The Ser529 site on the C-terminal transactivation domain of NFκB-p65 is often phosphorylated in response to the same stimuli that result in degradation of IκB. This phosphorylation improves the transcriptional activity of p65 but does not affect nuclear translocation or DNA binding. The K10-895.12.50 monoclonal antibody recognizes the phosphorylated serine (pS529) of human NF-κB p65 subunit.

Application: The metal-tagged antibody is designed and formulated for the application of Imaging Mass Cytometry™ (IMC™) using the Fluidigm Hyperion™ Imaging System on formalin-fixed, paraffin-embedded (FFPE) tissue sections.

Quality control: Each lot of conjugated antibody is quality control- tested by Imaging Mass Cytometry on tissue sections

Recommended concentration: For optimal performance it is recommended that the antibody be titrated for the desired application. Suggested initial dilution range: IMC-Paraffin: 1:25 to 1:100


References

Chang, Q. et al. "Staining of frozen and formalin-fixed, paraffin-embedded tissues with metal-labeled antibodies for Imaging Mass Cytometry analysis." *Current Protocols in Cytometry* 82 (2017): 12.47.1–12.47.8.

Giesen, C. et al. "Highly multiplexed imaging of tumor tissues with subcellular resolution by mass cytometry." *Nature Methods* 11 (2014): 417–22.

Safety

Use standard laboratory safety protocols. Read and understand the safety data sheets (SDSs) before handling chemicals. To obtain SDSs, go to fluidigm.com/sds and search for the SDS using either the product name or the part number.

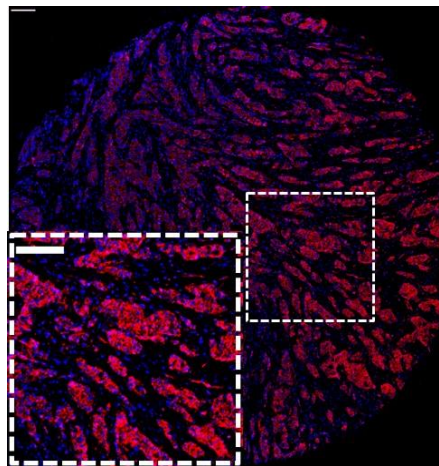

Human breast carcinoma (FFPE) stained with 166Er-anti-pNFkBp65 [S529] (K10-895.12.50) at a dilution of 1:50 (red pseudocolor) and iridium DNA intercalator (blue pseudocolor). Heat-mediated antigen retrieval was performed using Tris/EDTA buffer pH 9.

Scale bar size = 100 µm.

For technical support visit techsupport.fluidigm.com. | For general support visit fluidigm.com/support.

For Research Use Only. Not for use in diagnostic procedures.

Information in this publication is subject to change without notice. **Limited Use Label License:** The purchase of this Fluidigm Instrument and/or Consumable product conveys to the purchaser the limited, nontransferable right to use with only Fluidigm Consumables and/or Instruments respectively except as approved in writing by Fluidigm. **Patent and license information:** fluidigm.com/legal/notices. **Trademarks:** Fluidigm, the Fluidigm logo, Hyperion, Imaging Mass Cytometry, and IMC are trademarks and/or registered trademarks of Fluidigm Corporation in the United States and/or other countries. All other trademarks are the sole property of their respective owners. © 2020 Fluidigm Corporation. All rights reserved. 24-Apr-2020